
Light Therapy skin mask

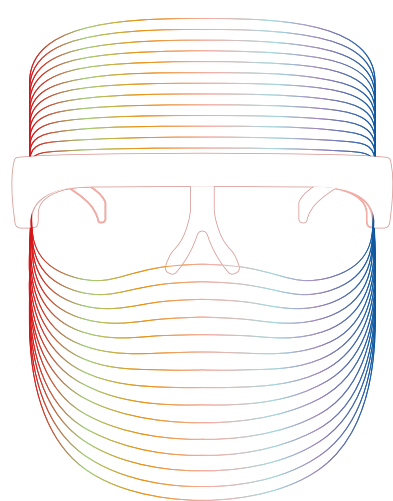User Manual

Please read the manual carefully to ensure safe and correct use.

Caution

Only use the mask wirelessly

- Use it on a clean skin, do not use harsh products such as acidic products, spot treatments and retinols before using the mask.
- Electric circuit inside the shield of the device should not be pressed or folded.
- Store the mask back in box to ensure the shield does not damage.

Principles of the product

This product uses 3 types of light sources including red, blue and orange.

Name of each component

1. USB charging cord
2. Nose piece
3. Switch (On, Off and changing the colours)

Product functions

The processes of LED optic therapy belongs to non-invasive painless skincare, and their effects are significant without harmful light. Hence, it is used widely by dermatologists. This product uses 3 LED light colours to improve your skin. Each colour assists different factors with different wavelengths:

Red light (620nm-630nm)	Assists blood circulation, accelerates metabolism, generates collagen, assists wrinkles and overall skin improvement.
Blue light (460nm - 470nm)	Helps improve acne and skin inflammation
Orange light (580nm-590nm)	Improves oxygen in the skin, accelerates the blood circulation and fades pigments. It is also useful to soften the skin.

Product specifications

Model	M series
Product name	Light Therapy skin mask
Product size	195*170*120mm
Product weight	70g
Material	ABS + PET + PVC
Colours	Red, Blue & Orange
Accessories	LED mask, manual and charging cable

How to use

LED Light Therapy can be used daily for optimal results

1. Plug the mask in using the charger provided
2. Cleanse the skin
3. Apply soothing and hydrating products if you wish (do not use harsh products such as acidic products, spot treatments and retinols before using the mask)
4. Push button for 3 seconds to TURN ON the shield
5. Press the button until desired LED light colour appears
6. Place the shield on your face for 10 to 20 minutes for optimal results (can be done daily)
7. Hold button down until shield POWERS OFF

Maintenance and storage

Maintenance:

1. Wipe the shield with a dry cloth after use, do not use water.
2. Always unplug after use and do not disassemble the product.

Storage:

1. If possible, store the mask in the original box to assist keeping the shape and not damaging the shield.
2. Store this product in a cool and dry place avoiding humidity and direct sun light.
3. Do not submerge the product in water or liquid.
4. Keep out of reach of children. Only suitable for adults.

Frequently asked questions (FAQs)

1. Can a minor use this product?

This product is only recommended for adults.


2. Does the light of the product have impact on the eyes?

Due to the unique design of this product, it illuminates light sources evenly on the face and is structurally designed to not harm the eyes. In the case of sensitivity, please use the product with your eyes closed.

3. Is there a limit to how often this product can be used?

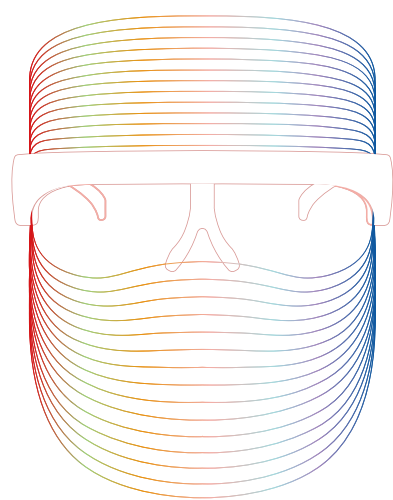This device is safe to be used everyday for 10-20 minutes per day.

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