LCD model (8 modes, 19 levels)

- It can be used for approximately 160 minutes after one hour of charging. The greater the intensity, the greater the power consumption will be. Please note that you can't use it during charging; it can only be used after charging.
- Place your feet on it before turning it on. Do not wear socks and remove the plastic bag. Do not use it while charging.
- The button on the main unit is like a button on the mat. Take off your socks and place your feet on the footprints. Press the "+" button to turn on the machine and press it twice to increase the intensity. The "-" button indicates reduced intensity or shutdown, which also indicates switching modes

What is EMS technology?

Electrical Muscle Stimulation (EMS), also known as Neuromuscular Electrical Stimulation (NMES), is the use of electrical impulses to induce muscle contractions. These pulses mimic action potentials in the central nervous system, causing muscle contractions. Studies have proven that EMS is a highly effective tool for positive muscle strengthening.

What is the purpose of EMS?

- Muscle strengthening
- Preventing muscle atrophy
- Muscle relaxation
- Relieving soreness
- Beauty conditioning
- Relaxing one's muscles

Troubleshooting

Conditions	Possible reasons	Solutions
Fail to connect power	The main unit battery is exhausted.	Please change the host computer.
	Used in a low or high-temperature environment.	It is recommended to use the ambient temperature between 15~35°C.
No electrical stimulation	The power is not connected.	Short press the on or off button until the LED light is on.
Current stimulation is too weak	Intensity setting is too low.	Press harder to increase intensity.
Current stimulation is too strong	Intensity setting is too high.	The compression is weakened and the intensity is reduced.
Current stimulation is suddenly interrupted	The remaining power of the main unit battery is exhausted.	Please charge the host computer.

Dangerous

Do not use this product with the following medical electronic instruments. Otherwise, it may cause accidents or malfunction of medical electronic instruments, and bring serious damage to the body. (Pacemakers; Defibrillators; Life-support equipment; Monitoring devices)

Note: Take off your socks when using the device, place both bare feet on the mat at the same time, and then press the Start button.


Do not use by the following individuals. Otherwise, it may result in accidents or skin and body problems.


Pregnant women, women who have just given birth, and individuals who are sensitive to electronic devices should not use this product. This includes: Those with a pacemaker or automatic implantable cardiac defibrillator (AICD). Individuals with heart disease. Those with high blood pressure.

EMS Foot Massager Instruction manual

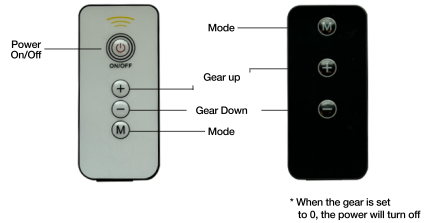
The packaging includes one of the remote controls shown above. Please note that the specific remote control provided may vary.

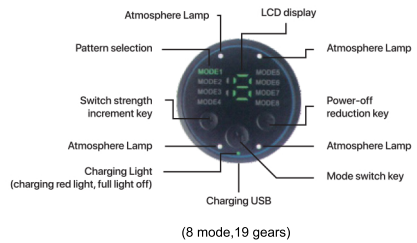Host usage instructions


- The stimulation intensity is divided into 19 levels (intensity 0~19 levels).
- When using it for the first time, please set a lower strength. After adapting to the low strength, you can adjust to other strengths.
- As the remaining battery capacity gradually decreases, the current stimulation will also gradually weaken. If you continue to use it, please charge the main unit promptly.

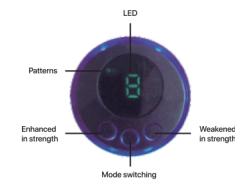(8 modes)

When the machine is switched to levels 1-19, it will automatically shut down every 15 minutes.

Charging: About 1-2 hours. The red light goes out when it is fully charged.

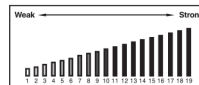Usage Tip: Please start with a lower intensity when using it for the first time, and adjust it after adapting to the low intensity.

Modes: There are 8 modes in total. For details, please refer to the following massager mode guide.

Technical Specification

Product name	EMS Foot Massager
Voltage	3.7V
Usage	A-F (8 modes) 19 levels
Frequency	1~100Hz
Working hours	15 minutes
Output	Maximum 35mA
Size	62x62x16
Weight	31
Battery capacity	150mAh

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